


AMM 215WTTP

Full HD 21.5-inch medical touchscreen

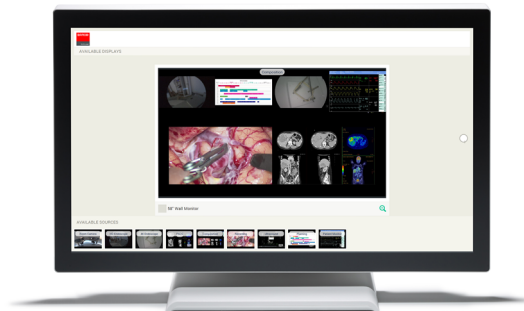

- **Accurate touchscreen**
- **Ergonomic and easy to operate**
- **Medically certified**

The AMM 215WTTP is a 21.5" touchscreen designed for medical professionals that demand accurate and reliable touch User Interface interaction.

Ergonomic interaction

The AMM 215WTTP has been built with ergonomics in mind. The 21.5" screen is easy to operate and the perfect size and format for a variety of applications, including the digital OR and medical carts. Thanks to its Full HD resolution, clinicians have access to rich and detailed information.

The display features Projected Capacitive Touch screen technology with multi- and single-touch functionality. It allows clinicians to smoothly interact with the user interface on the monitor, even when wearing surgical gloves.

Designed for medical environments

The AMM 215WTTP is medically certified and fit for use in a highly demanding, mission-critical environment, such as the operating room. The bezel-free design allows for easy cleaning and disinfection, a major prerequisite in all surgical and clinical environments.

Easy to install

The AMM 215WTTP allows for easy installation onto a variety of computing platforms. It is Microsoft® Windows® HID (Human Interface Device) compatible if you use the USB touchscreen interface. No additional software driver is required for general operation of the touchscreen.

- Project Capacitive touchscreen with multi-touch and gesture input application

- Full High Definition resolution
- Medically certified
- Bezel-free design
- HID complaint for easy installation & integration on computing devices
- Easy cleaning
- Table stand allowing for wide tilt angle
- VESA mounting holes on touch monitor unit for mounting on surgical boom arms, monitor arm and medical carts
- Mounting holes on the table stand for various mounting options
- Multilingual OSD user control

PRODUCT SPECIFICATIONS**AMM 215WTTP**

| | |
|-------------------------------|---|
| Screen technology | TFT AM LCD / LED Backlight |
| Active screen size (diagonal) | 21.46" inches (545.2 mm) |
| Active screen size (H x V) | 476 x 267 mm (18.7 x 10.5") |
| Aspect ratio | 16:9 |
| Resolution | 1920 x 1080 |
| Pixel pitch | 0.248 x 0.248 mm |
| Color support | 16.7 million colors |
| Viewing angle | 178° Hor / 178° Ver |
| Luminance | 250 cd/m ² (typ) |
| Contrast ratio | 1000:1 (Typical) |
| Touch sensor | Projected Capacitive Touch Technology (PCT) 10 points |
| Video input signals | 1xDVI-D, 1xVGA (D-Sub 15pin) |
| Video output signals | N/A |
| Video formats | Up to 1920x1080, 60 Hz 31.47-82.3 kHz (Horizontal) 56-75Hz (Vertical) |
| Interface | USB 2.0, type B |
| Touchscreen driver | Windows, Embedded Linux, Apple Mac OS |
| External power supply | AC input: 100 – 240 VAC / 50-60 Hz DC output: +12VDC / 4.2 A |
| Power consumption | Operating: 27W Sleep: 2 W -Off : 1 W |
| Dimensions (W x H x D) | 518 x 314 x 46 mm 518 x 352 x 219 mm (incl. stand) |
| Net weight display | 7 kg (approx.) |
| Net weight packaged | 10 kg (approx.) |
| Mounting standard | Removable display stand 100x100 VESA mount |
| Temperature | Operating: 0° -40° C Storage: -20° -60° C |
| Humidity | Operating: 20% – 80% RH, non-condensing Storage: 10% – 80% RH, non-condensing |
| Regulation Compliance | CE (EN60601-1, 3rd edition) and EMC (EN60601-1-2) UL/cUL (ANSI/AAMI ES60601-1, CAN/CSA-C22.2 No. 60601-1) FCC (Part 15 Class B) CCC (GB4943.1, GB9254, GB17625.1) PSE |
| Green compliance | RoHS, REACH |
| Warranty | 18 months |

Last updated: 29 Nov 2021

Technical specifications are subject to change without prior notice. Please check www.barco.com for the latest information.